



Insertion of TruSteel™ Infusion Set

This guide is a summary only and intended as a quick reference guide only. Please refer to package insert for full instructions and precautions.

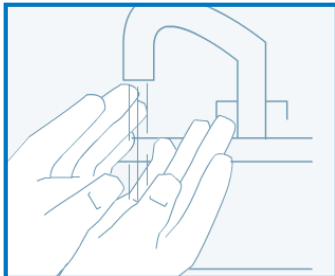


Wash and dry your hands thoroughly.

Thoroughly clean the insertion site and allow to dry. Use of **Adhesive Wipes** will clean the site and help the set to stay in place.

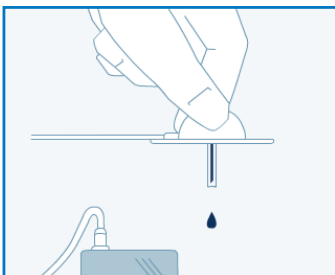


Connect the infusion set line to the tubing connector of the cartridge and then connect the needle part of the infusion set to the tubing.

Select Fill Tubing to push the insulin through the line until you maintain a steady flow of insulin for at least 10u. Ensure there are no air bubbles.

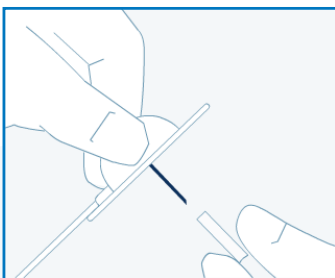


Remove the adhesive tape backing and then the needle protector from the set.

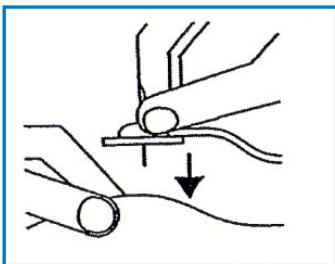

Hold the infusion set just above the needle with one hand. Pinch the skin at the injection site with the other hand.

Insert the infusion set to its full length at a 90° angle (straight in).

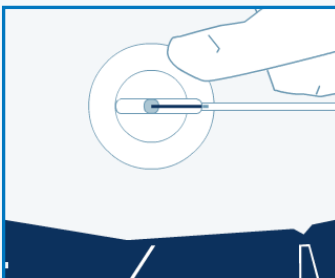

Ensure the material is properly placed and adequately covers the injection site. Remove the adhesive back taping from the connection point and secure firmly onto the skin.

There is no cannula fill required to be entered into the pump for this infusion set.

Your Choice. Your Life. Your Freedom.®

24/7 Technical Support 0508 634 103  
0800 500 226 | nzms@nzms.co.nz | nzmsdiabetes.co.nz

