


# Insertion of AutoSoft™ 30 Infusion Set

This guide is a summary only and intended as a quick reference guide only. Please refer to package insert for full instructions and precautions.

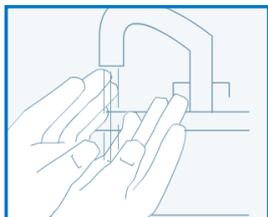


Wash and dry your hands thoroughly.


Thoroughly clean the insertion site and allow to dry. Use of [Adhesive Wipes](#) will clean the site and help the set to stay in place.




Pull the red tape to remove the seal.

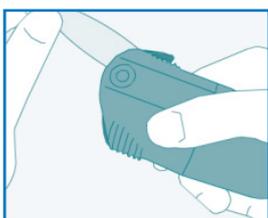

Remove the sterile paper from the bottom of the set housing.

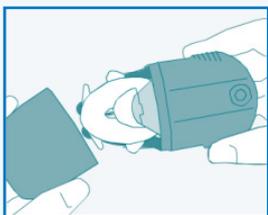

Hold the lined indentations with one hand pull the lid off with the other.

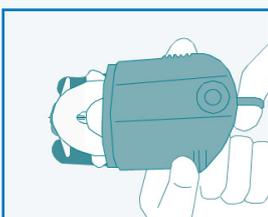


With one hand, place your fingers on the indentations on either side of the housing.

With the other hand pull back the spring softly until you hear a click. The needle protector will automatically lift up or peel off.




Carefully bend back the needle protector. This will also remove part of the adhesive tape. Ensure that the tape is not stuck to the needle.

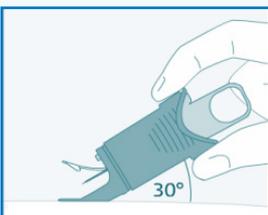

Place your index finger on the top release button and your thumb on the bottom release button. Position AutoSoft™ 30 with its legs level (flat) on the skin to ensure an insertion angle of 30°.




Gently squeeze on the top and underside of the housing to insert the set.

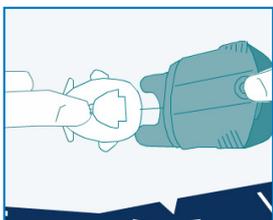


Secure the infusion set with one finger in front of the clear window.

Remove the insertion device and introducer needle by pulling gently straight back.

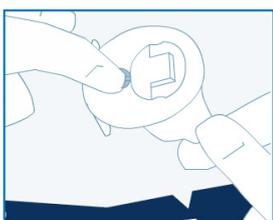


Place one finger in front of the clear window. Gently remove the backing tape from the other half of the site. Hold skin tight and gently press outside of adhesive tape to secure edges.

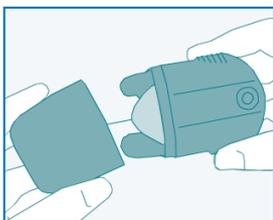

Replace the lid back on for safe disposal.



Connect the line to the cartridge in the pump and fill tubing. Ensure there are no air bubbles.

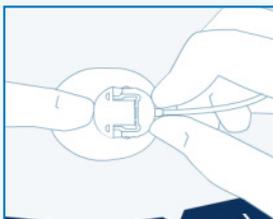

Connect the line to the site and make sure an audible click is heard and fill the cannula.

**REMEMBER TO FILL NEW CANNULA WITH 0.70 UNITS**

Your Choice. Your Life. Your Freedom.®

24/7 Technical Support 0508 634 103    
0800 500 226 | [nzms@nzms.co.nz](mailto:nzms@nzms.co.nz) | [nzmsdiabetes.co.nz](http://nzmsdiabetes.co.nz)

NZDI\_TAN\_017

